

LABORATORY. Our range of laboratory refrigerators, freezers and fridge-freezers offer the safe and reliable storage of volatile liquids, reagents, chemicals, poisons and samples, while our ULT freezers and laboratory incubators have been developed for the safe storage of sensitive specimens such as those for biomedical research. In an exclusive deal, we also offer a suite of ice-makers from Lec Medical by Hoshizaki to compliment our laboratory range.



IST147

Safety

Safety inside a laboratory is fundamental so when choosing a laboratory fridge or freezer spark-free is the safest option. By removing the fan and internal light this eliminates the risk of explosion within the interior.

Choice

A wide range of sizes and variants in all laboratory products provides a great deal of choice, with many fridges and freezers having matching options available.

Control

Easy to use external control panels and digital temperature displays ensure continuous monitoring without having to open the door and disturb the internal atmosphere.

Security

All Lec Medical laboratory refrigeration products are fully lockable preventing unauthorised access at all times. Our incubator range includes a password protect function preventing any unauthorised adjustment to the control and alarm thresholds.

Environment

Low energy consumption and being CFC free makes all Lec Medical laboratory appliances kinder to the environment.



LABORATORY FRIDGES. The laboratory fridges from Lec Medical are available in five different sizes. From 82 Litres to 475 Litres with the four smaller options having matching laboratory freezers available.

Control and Efficiency

The integral control panel on the outside of the fridge sits flush with the door and its easy to use design makes checking the internal temperature straightforward. As there is no need to open the door to check the temperature, it keeps the internal air temperature more stable.

Automatic defrost on all of the fridges ensures more efficient cooling.

Security

The spark-free interior reduces the risk of internal explosion.

Security is essential in laboratories so all Lec Medical fridges come complete with door locks and two keys ensuring that only authorised individuals are able to access them.

A temperature alarm, which sounds if the temperature deviates outside the set parameters, comes as standard on all laboratory fridges reducing the chance of ruining the contents of the refrigerator.

Flexibility

Reversible doors* give greater flexibility where either space is at a premium or positioning the product is specific.

Adjustable feet are available on most models making it easier when positioning the product on uneven surfaces.

All of the laboratory fridges come with removable shelves giving maximum storage capacity.

Environmentally Friendly

All our products are CFC free with many being HFC free too. Lec Medical refrigeration products are some of the most efficient products in the market using both efficient refrigerants and high specification insulation to help drive down the energy cost for our customers.

LR207

SKU	444448568
Product dimensions (mm)	502w x 540d x 660h
Packed dimensions (mm)	520w x 595d x 675h
Product weight (kg)	27.5
Packed weight (kg)	30.5
Net capacity (litres)	67
Gross capacity (litres)	82

LR307

SKU	444449326
Product dimensions (mm)	495w x 580d x 830h
Packed dimensions (mm)	560w x 655d x 850h
Product weight (kg)	31.5
Packed weight (kg)	35
Net capacity (litres)	103
Gross capacity (litres)	107

LR507

SKU	444448514
Product dimensions (mm)	558w x 600d x 863h
Packed dimensions (mm)	600w x 610d x 900h
Product weight (kg)	34
Packed weight (kg)	35.5
Net capacity (litres)	145
Gross capacity (litres)	153

LR907

SKU	444448516
Product dimensions (mm)	595w x 600d x 1550h
Packed dimensions (mm)	625w x 700d x 1610h
Product weight (kg)	81.5
Packed weight (kg)	84
Net capacity (litres)	285
Gross capacity (litres)	286

LR1607

SKU	444448569
Product dimensions (mm)	685w x 910d x 1688h
Packed dimensions (mm)	770w x 950d x 1810h
Product weight (kg)	130
Packed weight (kg)	140
Net capacity (litres)	441
Gross capacity (litres)	475



LR207




LR307



LR507



LR907

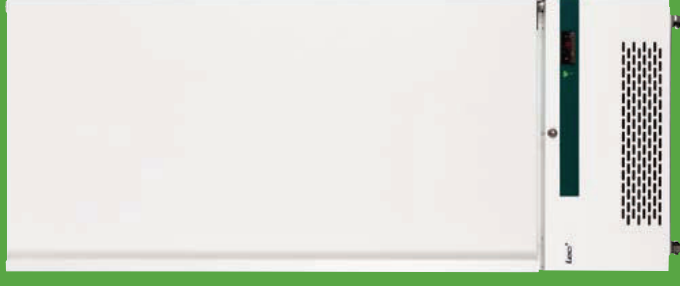

LR1607